molift[®]

EASY

Molift Easy Ambulating Vest

Molift offers a broad range of slings for all movement needs. The slings are designed to distribute weight and pressure as comfortably as possible, so the slings vary both in model and size. Molift focuses on safety and comfort, and recommends individual fitting of the slings. Colour codes help you to find the correct size. We design and produce our slings in Norway, in accordance with the same quality constraints as the Molift lifts, and all slings are tested and approved in accordance with NS-EN ISO 10535. Molift recommends using Molift slings with Molift lifts.

Molift Easy Ambulating Vest


Size	Colour code	Article no.
S	Red	Art. no 3036100
M	Yellow	Art. no 3036200
L	Green	Art. no 3036300

- a sling from ***molift***[®]
designed for life


USER MANUAL

Molift Easy Ambulating Vest

Ambulating Vest has been specially developed for rehabilitation, and is intended for use with Molift's 4 point suspension system. It is also very suitable for toilet visits. The high belt is extra padded so that pressure and load are evenly distributed. The belt is equipped with double belt buckles for diagonal tightening, so that you achieve optimum support and stability. The groin straps are well padded and prevent the Ambulating Vest sliding up. Self-tightening suspension loops prevent the user sliding or falling if his or her legs give way. This system also allows the lift to be raised or lowered to adjust the load on the user's legs. All of these details permit a high degree of comfort and safety for the user.

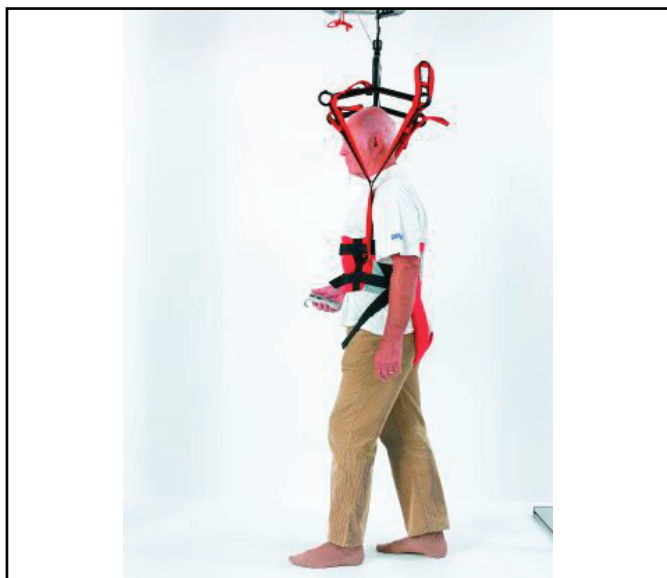The buckles in the front allow the user to use the Ambulating Vest alone.

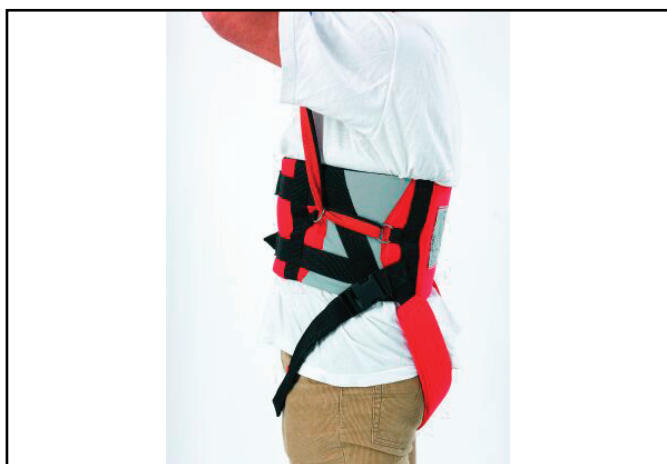Molift Ambulating Vest is made from grey and red polyamide material (Cordura), with polyester tapes, padded with polyester foam, and can be washed.

Torn, worn or in any way damaged slings/straps can break, and lead to injury. Only use slings that are in good condition. Destroy and discard old, useless slings.

*Lifting and moving patients always involves a certain degree of risk. **Read the lift's user manual carefully!** It is important to understand and be familiar with the equipment, and only trained personnel should use it. Ensure that the sling fits the lift you have chosen, and that it is correct for the situation it is to resolve. Show caution and care when you use the lift and its accessories because you, as a carer, are responsible for the safety of the patient.*

Molift Easy Ambulating Vest is designed to fit the Molift 4-point suspension and is intended as an accessory to Molift lifts. The lift's user manual indicates which slings can be used with it.

Using the Molift Easy Ambulating Vest

Molift Easy Ambulating Vest is usually used in the following situations:

- Standing, walking and balance training with the Molift ceiling-mounted lift or Molift Partner 3 (see accessories, user manual)
- Toilet visits using the Molift ceiling-mounted lift

Fitting the sling

- Lean the user forward while you hold one arm round him or her. Lead the seat down the user's back to waist height.

Remember the correct working position to protect your back!


- Fix the waist belt and tighten.

i Remember that the upper body is slimmer when a person is standing erect, so it is always important that the buckles are adjusted after standing for maximum safety and stability.

- Position the groin strap under the thigh, and fix the rear buckle on the side of the back.

- Molift Easy Ambulating Vest is flexible and is very adjustable. Always check that you find the best solution for your user.
- The lift suspension points must always be arranged across the user as shown to the left in the illustration below.

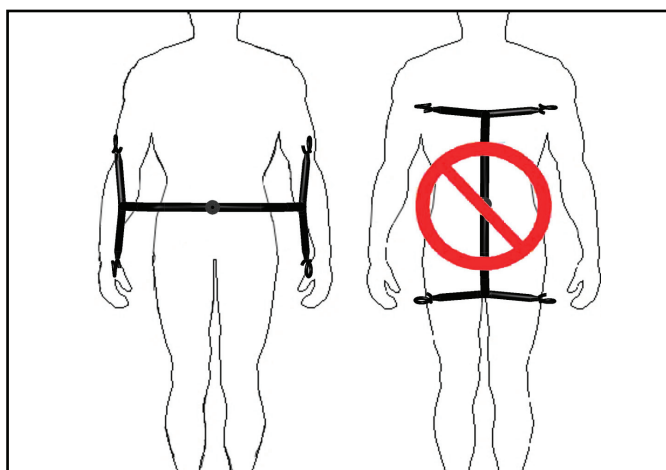- Hook the loops onto the lift hooks.

! Ensure that the suspension loops are affixed correctly to the suspension points so that there is equal load on the shoulder and leg loops.

- If you are lifting from a wheelchair, apply the wheelchair brake so it does not move. Lift, at the same time as the user raises him- or herself erect. Ensure that you do not lift too high, so that the user hangs in the straps, but feels that he or she is standing safely and comfortably.
- Lower the lift as the user sits carefully down into the chair, take the loops off the hooks and remove the lift. If moving to a wheelchair, the brakes must be on to prevent the wheelchair from moving.

! Remember the correct working position to protect your back!

- Remove the groin strap and remove it from below the user's thigh. Remove the belt support the user with one hand. Pull the sling up from behind the user's back. Lean the user back in the chair.


! If you pull the sling up hard, the user may fall forwards and be injured. Always hold one hand around him or her when removing the sling.

Washing instructions

Molift Easy Ambulating Vest is made from grey and red polyamide material (Cordura), with polyester tapes, padded with polyester foam, and can be washed at 85° C and autoclaved at 85° C for 30 min.

Size guide for sling sizes

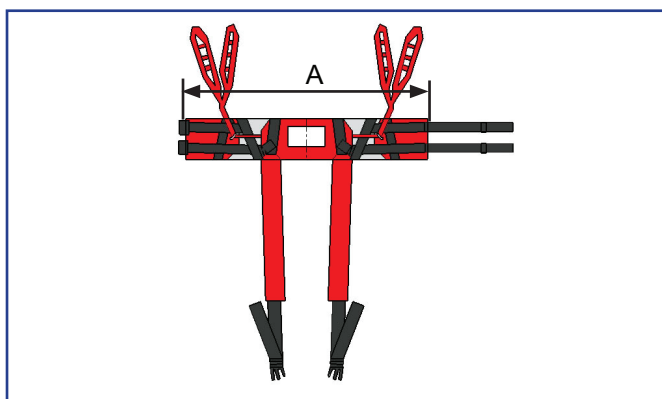Size	Colour code	Recommended user weight (Kg)	Waist size A (cm)
S	Red	25-50	60-75
M	Yellow	45-95	75-100
L	Green	90-160	100-150

The choice of sling size is determined by the user's weight and body shape/size.

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